

Weight limit 265 lbs

CD101XT Summit™ Extra Tab




Customer Service
(800) 819-5980 or (208) 429-0026
www.coyotedesign.com



Coyote Design & Mfg., Inc.
7450 W Lemhi Street, Boise, ID 83709

MADE IN USA

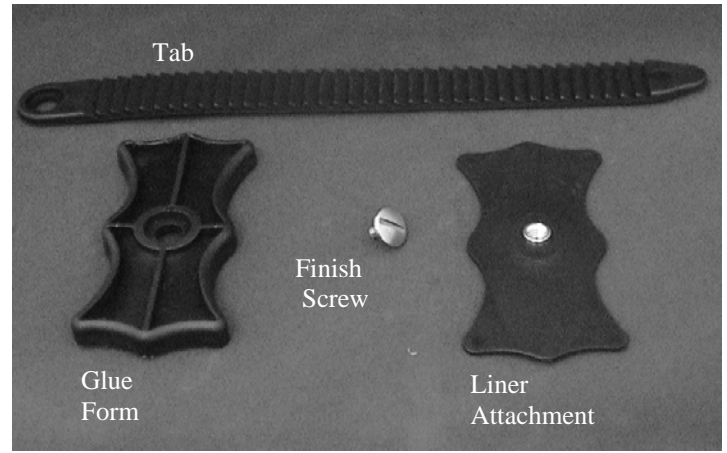

Gluing Liner Attachment to Liner

Correct fabrication of the Summit™ Lock requires use of Coyote® Quick Adhesive. The Quick Adhesive sets extremely fast. If not glued fast enough it may not properly adhere to the liner. The glue can be refrigerated to increase the work time. When using, work quickly.

The Quick Adhesive is sold separately.

CD4120—20cc cartridge

CD4180—180cc Cartridge


* CD4180 requires a dispensing gun and mix tips



Patient steps into socket. Liner is marked at top edge of hole in socket.



Roll liner over a piece of plastic to pre-stretch the liner and to create a flat working surface.

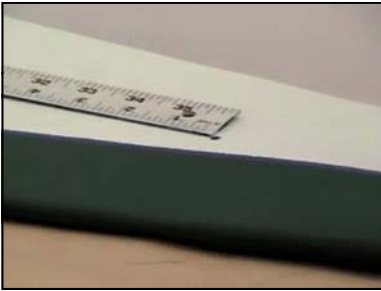

Measure 1/2" to 1" above mark on liner. This mark represents the bottom of the Bowtie attachment.



Attach liner attachment to the glue form with the finish screw.



Clean the underside of the liner attachment with acetone.



Mark around the liner attachment to ensure exact placement.



Run a bead of Quick Adhesive around the edge of the liner attachment. Smear the bead into a thin layer on the attachment.



Quickly set liner attachment in marked location.



Place small weight on top of the glue form to hold in place. Allow to set for 15 minutes.



Remove the glue form and secure the tab to the liner attachment with the finish screw. Use Loctite® 242 to ensure permanent attachment.